

SOMATEX®
Medical Technologies GmbH

Rheinstraße 7d D-14513 Teltow Germany

Tel. +49 (0) 3328 30 76-0 Fax +49 (0) 3328 30 76-99 info@somatex.com

DUO System PREMIUM – for easy, single-handed operation

The DUO System PREMIUM from SOMATEX® is a correctable localization system for mammography. It offers doctors and patients top quality, absolute reliability and highest comfort for the preoperative marking of suspicious tissue.


DUO System PREMIUM is particularly userfriendly for ultrasound treatments as it can be used single-handedly. This allows the user to make markings simultaneously to guiding the ultrasound transducer.

The removable handle of the DUO System PREMIUM constitutes a real contribution to the mammography procedure: By means of punctures through a breadboard, it allows for a new positioning of the compression plates without having to remove the cannula.

If a repositioning is required, the marker wire of the DUO System PREMIUM can be easily retracted into the cannula thanks to its special form. The wire can then be released again in the intended and correct position. The double arches provide firm anchorage in the tissue and therewith prevent dislocations in any direction. The varying features of the wire allow to recognize the depth position of the marked lesion.

DUO System PREMIUM offers easy and precise handling as well as a gentle application and ensures optimum visibility for the ultrasound procedure.



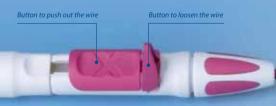



Removable cannula for puncturing through breadboard

Varying threading of the wires every 50 mm allows for depth orientation.

Advantages Product options

Use


DUO System PREMIUM serves for the preoperative marking of nonpalpable, suspicious mammary lesions and simplifies intraoperative localization.

Advantages at a glance

- Single-handed operation
- Removable cannula for puncturing through breadboard
- Correctable marker wire
- Firm anchorage in the tissue by means of DUO arches
- Sharp wire tips for hard lesions
- High flexibility of the marker wire
- Depth orientation at the cannula and wire
- Ultrasound marking and ultrasound cut at the cannula tip

with cm marking

Fixation piece allows for the safe preoperative fixation of the wire

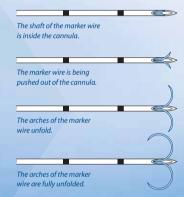

Correctable arches

Ultrasound marking

Ultra.

Ultrasound cut

DUO System PREMIUM			
Order no.	Gauge	Diameter	Length
271 770	19.5	0.95 mm	50 mm
271 771	19.5	0.95 mm	90 mm
271 773	19.5	0.95 mm	120 mm

DUO System PREMIUM