

Tuloc[®]Premium

Correctable localization system for easy,
single-handed operation



SOMATEX[®]
MEDICAL TECHNOLOGIES GMBH

www.somatex.com

Tuloc® Premium

for easy, single-handed operation

Tuloc® PREMIUM is a correctable localization system. It is intended for the preoperative marking of non-palpable, suspicious breast lesions and facilitates intra-operative localization for the surgeon. Tuloc® PREMIUM offers doctors and patients top quality, absolute reliability and highest comfort.

Tuloc® PREMIUM is particularly user friendly for ultrasound treatments as it can be used single-handedly. This allows the user to make markings simultaneously to guiding the ultrasound transducer. The removable handle of the Tuloc® PREMIUM constitutes a real contribution to the mammography procedure: By means of punctures through a breadboard, it allows for a new positioning of the compression plates without having to remove the cannula.

If a repositioning is required, the marker wire of the TULOC® PREMIUM can be easily retracted into the cannula thanks to its special form. The wire can then be released again in the intended and correct position. The double arches provide firm anchorage in the tissue and there-with prevent dislocations in any direction. The wire of the Tuloc® has markings for secure location for the user. Two proximal markings indicate whether the arches of the wire are inside or outside of the cannula. The distal markings at distances of 2 cm give information on the distance to the tumor.

Tuloc® PREMIUM offers easy and precise handling as well as a gentle application and ensures optimum visibility for the ultrasound procedure.

Use

Advantages


Product options

Advantages at a glance

- Single-handed operation
- Removable cannula for puncturing through breadboard
- Correctable localization wire
- Monofilament wire provides for best possible form and pressure stability
- Distal depth markings show information on the distance to the tumor
- Extreme sharpness for precise and atraumatic puncture
- Special material and shape of double loop secure good visibility for ultrasound
- Ultrasound marking and ultrasound cut at the cannula tip
- Fixation piece allows for the safe preoperative fixation of the wire

Use


- Tuloc® PREMIUM is intended for the preoperative marking of non-palpable, suspicious breast lesions and facilitates intraoperative localization for the surgeon.

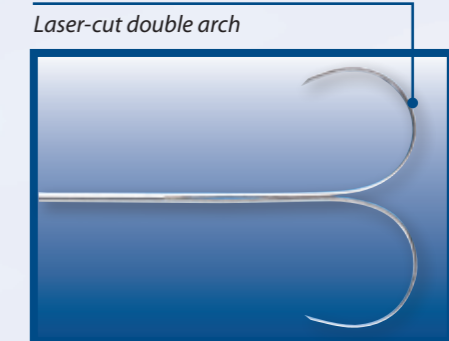

This marking indicates that the arches are still inside the cannula

This marking shows that the arches are completely unfolded

Removable cannula for puncturing through breadboard

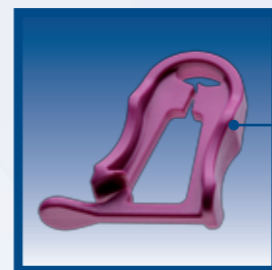

Fixation piece

4 cm to the distal tip of the wire

2 cm to the distal tip of the wire

Markings every 2 cm provide information on the distance to the distal tip of the wire

Ultrasound cut

Correctable arches

Ultrasound marking

Cannula with cm marking

Slider to put out the wire

Lever for loosening the wire

Tuloc® PREMIUM

Order No.	Gauge	Diameter	Length
271 880	20 G	0,95 mm	50 mm
271 881	20 G	0,95 mm	90 mm
271 883	20 G	0,95 mm	120 mm

SOMATEX[®]

MEDICAL TECHNOLOGIES GMBH

Rheinstr. 7d
D-14513 Teltow
Germany

Tel. +49 (0) 3328 30 76-0
Fax +49 (0) 3328 30 76-99
info@somatex.com
www.somatex.com