


Tumark® Q

Clip marking system

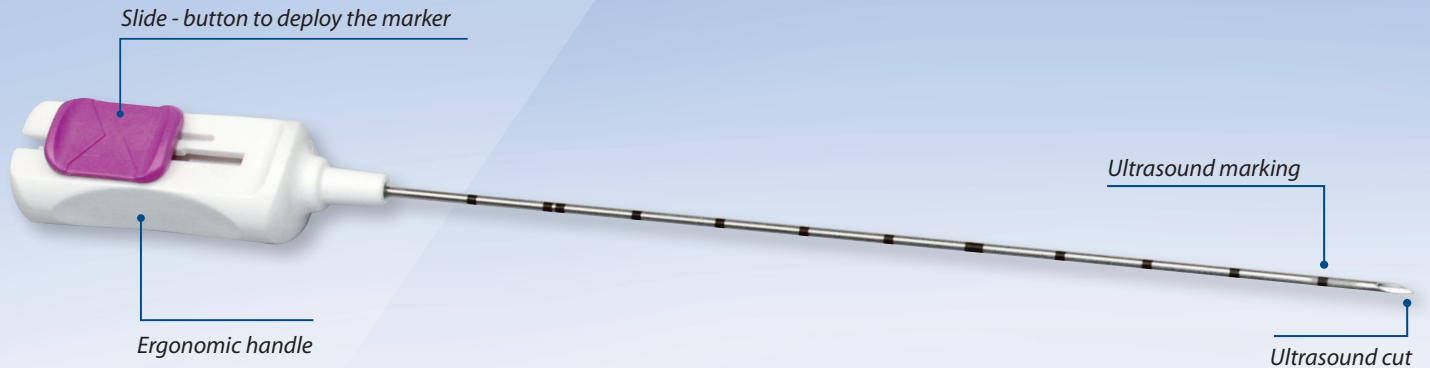

Product description:

The Tumark® Q serves for percutaneous marking in the pectoral tissue. Among other things, the marker is implanted in lesions prior to or during chemotherapy and the location of a biopsy sampling point. The marker is visible under mammography, ultrasound and MRI*.

* Only the marker of Tumark® Q is MRI compatible.


Advantages at a glance:

- Q-merged form enables firm anchorage in the breast tissue and increases visibility in all positions
- Visibility of the marker in ultrasound, mammography and MRI
- Marker is MRI compatible up to 3 Tesla
- Marker is approved as implantable material (Nitinol)
- Marker is supplied preloaded
- Ergonomic handle for one single-hand operation
- Simultaneous use with Tumark Professional (U-Marker) enables better differentiation during mammography

The marker is supplied preloaded



Clip-marker Q-shape

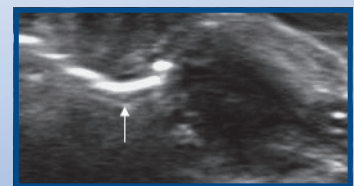

Ultrasound



Mammography

Article Number	Gauge	Diameter	Length
271 500	18 G	1,2 mm	120 mm

Sales Office:
 Hohenzollerndamm 150/151 Haus 2 · D-14199 Berlin · Germany
 Tel.: + 49 (0) 30 319 82 25-00 · Fax: + 49 (0) 30 319 82 25-99
 info@somatex.com · www.somatex.com

SOMATEX®
 MEDICAL TECHNOLOGIES GMBH

999953V1